

Instruction manual and data sheet bPCA-3000-05-10-800-x

Photoconductive THz antenna for laser excitation wavelengths $\lambda \sim 800$ nm

PCA – Photo Conductive Antenna


Table of contents:

1. Spectral performance.....	2
2. Antenna parameters.....	3
3. Antenna design.....	4
4. Order information.....	5




1. Spectral Performance

763 bPCA 3000-05-10 #05

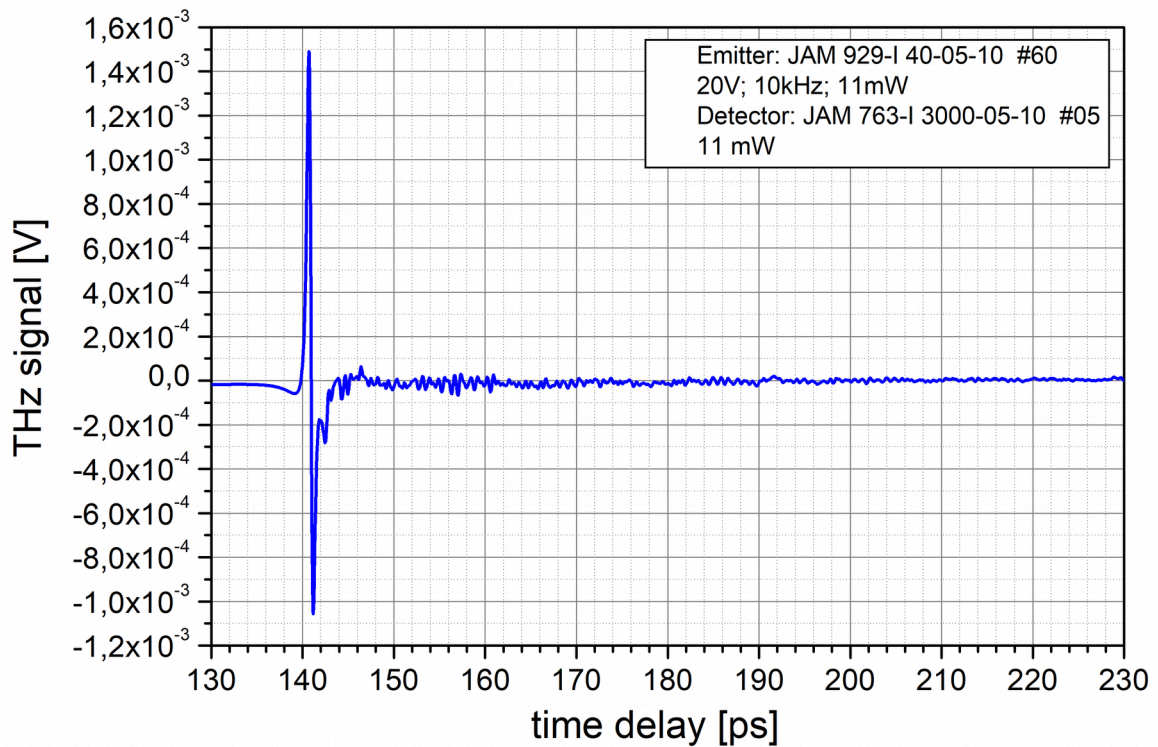

763 bPCA 3000-05-10 #05

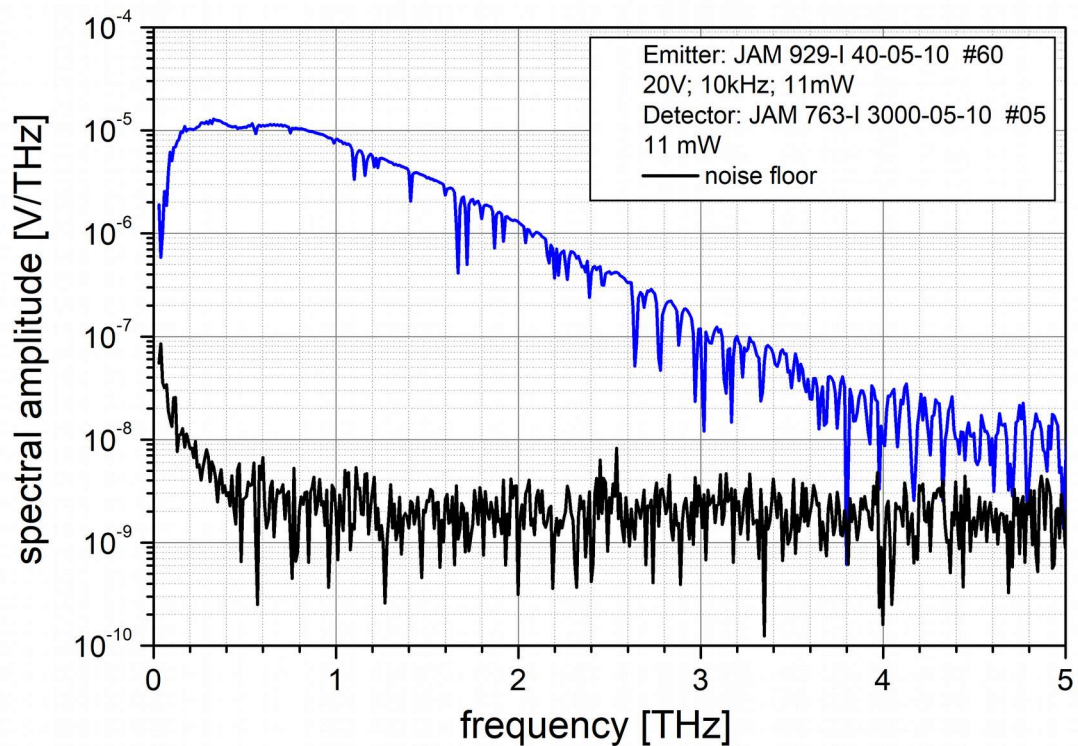

2. Antenna parameters

Parameter	minimum ratings	standard	maximum ratings
Dark resistance	500 MΩ	1 GΩ	2 GΩ
Voltage		20 V	25 V
Optical mean power @ 50 – 100 MHz repetition rate		10 mW	20 mW
Pulse fluence		250 μJ/cm ²	500 μJ/cm ²

Attention: The F-number of the optical lens focusing the laser beam onto the antenna gap must be larger than a certain value to avoid too high pulse fluency. This means, that the minimum diameter of the focused beam waist must be about 120 % of the gap distance g . For a Gaussian beam the minimum focus length f_{\min} of the optical lens can be estimated as

$$f_{\min} = \frac{0.3 \cdot \pi \cdot g \cdot D}{\lambda}$$


with g – gap distance of the antenna

λ - laser wavelength


D – diameter of the laser beam hitting the focusing lens.

For $\lambda = 0.8 \mu\text{m}$ and $g = 5 \mu\text{m}$ the minimum possible F-number of the lens is $f_{\min}/D = 1.9\pi$.


3. Antenna design

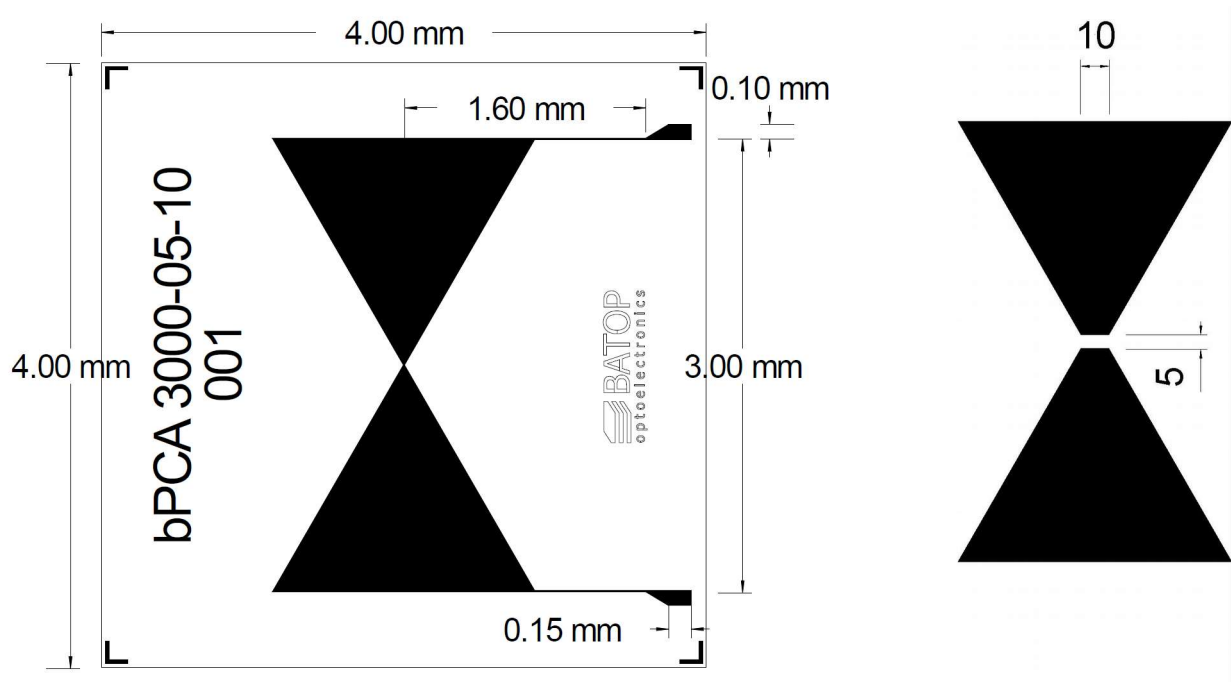

antenna dimensions in μm

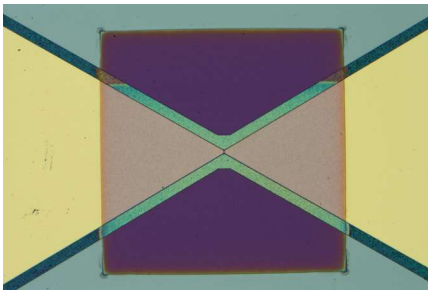


bPCA 3000-05-10-800

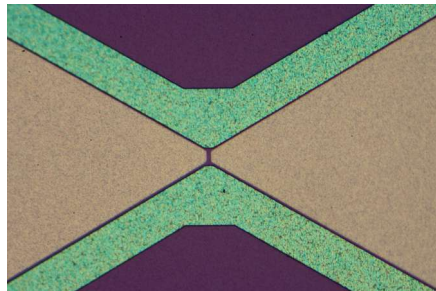

bPCA 3000-05-10-800 (detail)

4. Order information

bPCA-3000-05-10-800-x Photoconductive antenna
 length $l = 3000 \mu\text{m}$
 gap $g = 5 \mu\text{m}$
 width $w = 10 \mu\text{m}$
 laser wavelength $\lambda = 800 \text{ nm}$

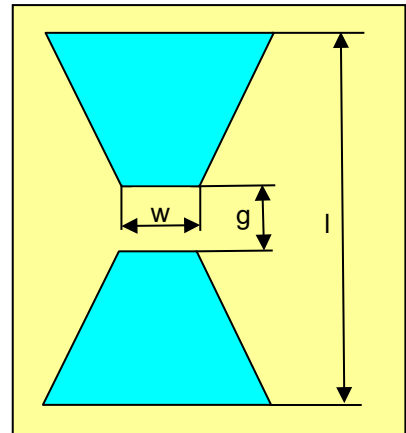

x denotes the type of mounting as follows:

- x = 0** unmounted chip 4 mm x 4 mm with 2 bond contact pads
- x = h** mounted on an Al disc with 25.4 mm \varnothing and [hyperhemispherical silicon substrate lens](#), 1m coaxial cable with BNC or SMA connector
- x = a** mounted on an Al disc with 25.4 mm \varnothing and [aspheric focusing silicon substrate lens](#), 1m coaxial cable with BNC or SMA connector
- x = c** mounted on an Al disc with 25.4 mm \varnothing and aspheric collimating silicon substrate lens CL-12 for 12 mm THz beam diameter, 1m coaxial cable with BNC or SMA connector
- x = c-f** [fiber coupled antenna](#) with collimating silicon substrate lens
- x = l** with [aspheric focusing optical lens](#) for free space laser excitation

For information about THz beam guiding possibilities please [click here](#)